



# **UE8** UPPER BODY ERGOMETER

**REDEFINING ACCESSIBILITY WITH THE INNOVATIVE SWIVEL DRIVE TRAIN**



**MOVEMENT  
FOR  
EVERY  
BODY™**

# INCLUSIVE FOR ALL USERS.

Enhance your lifestyle with the NuStep UE8 to provide an upper body cardio and strength workout – on one machine.

Low-inertia start up.

Bi-directional arm movement with 15 levels of resistance.

Intuitive swivel arm adjustment provides for a seated, standing, or from a wheelchair workout.

Rear power access – two power access points for placement flexibility.



Dual touch screen console with HD display and a user friendly interface.

Front power access – two power access points for placement flexibility.

Integrated sensor technology to capture power and symmetry data.

Low step through design for easy access.

Adjustable 360° swivel seat with reclining capabilities.

# ACCESSIBLE WORKOUT EXPERIENCE.

The UE8's unique swivel drive train design allows direct access from three positions – seated, standing and from a wheelchair – without requiring seat removal.




**Seated**  
The low step through and 360° swivel seat promotes ease of access for seated workouts.

**Standing**  
The swivel drive train easily adjusts up, down, right and left to accommodate users who want a standing workout.



## GET ON AND GO.

**Intuitive Touchscreen Console.** The UE8 console is designed to be low maintenance, and simple to use for every user.

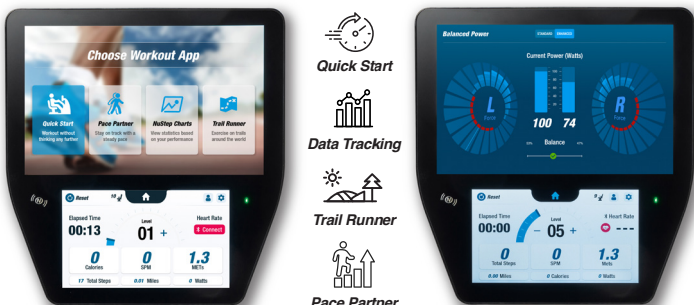

**Plus!**  
**Balanced Power** – Integrated sensor technology measures upper body power and symmetry in real time.



**Wheelchair**  
Direct access from a wheelchair is also made possible from the swivel drive train – without removal of the seat – providing an independent workout.

## SPECS

**Machine Footprint:**  
• Straight Only 74" x 32"  
• Right or Left Use Configuration: 74" x 54"  
• Both Right and Left Configuration: 74" x 76"

**Unit Weight:**  
• UE8PRO: 310 lbs.  
• UE8MAX: 320 lbs.

**User Weight:**  
• UE8PRO accommodates users up to 500 lbs.  
• UE8MAX accommodates users up to 600 lbs.

**Fits Users:**  
From 4'6" to 6'7" in height.

**Power:**  
Medical grade 24V DC power supply included.

**Warranty** The UE8PRO and UE8MAX are covered under a limited 10-year frame, 3-year parts and 1-year labor warranty (US). Product specifications are subject to change.



**MOVEMENT FOR EVERY BODY™** (800)322.2209 | (734) 769.3939 | [nustep.com](http://nustep.com)

Contact sales for more information about our accessories equipment! [sales@nustep.com](mailto:sales@nustep.com)